



# **International Rules for Seed Testing 2025**

**Validated Seed Health Testing Methods**

**7-002a: Detection of *Alternaria radicina* in *Daucus carota* (carrot) seed by blotter method**

**Including changes and editorial corrections adopted at the  
Ordinary General Meeting 2024 in Cambridge, United Kingdom**

**Effective from 1 January 2025**

## **Validation reports**

See References. Copies are available by e-mail from the ISTA Secretariat at [ista.office@ista.ch](mailto:ista.office@ista.ch).

Please send comments, suggestions or reports of problems relating to this method to the ISTA Seed Health Committee, c/o ISTA Secretariat.

## **Disclaimer**

Whilst ISTA has taken care to ensure the accuracy of the methods and information described in this method description, ISTA shall not be liable for any loss or damage, etc. resulting from the use of this method.

## **Safety precautions**

Ensure you are familiar with hazard data and take appropriate safety precautions, especially during weighing out of ingredients. It is assumed that persons carrying out this test are in a laboratory suitable for carrying out microbiological procedures and familiar with the principles of Good Laboratory Practice, Good Microbiological Practice, and aseptic techniques. Dispose of all waste materials in an appropriate way (e.g. autoclaving, disinfection) and in accordance with local health, environmental and safety regulations.

## **Note on the use of the translations**

The electronic version of the International Rules for Seed Testing includes the English, French, German and Spanish versions. If there are any questions on interpretation of the ISTA Rules, the English version is the definitive version.

Published by  
The International Seed Testing Association (ISTA)  
Richtiarkade 18, CH-8304 Wallisellen, Switzerland

© 2025 International Seed Testing Association (ISTA)

Online ISSN 2310-3655

All rights reserved. No part of this publication may be reproduced, stored in any retrieval system or transmitted in any form or by any means, electronic, mechanical, photocopying, recording or otherwise, without prior permission in writing from ISTA.

## 7-002a: Detection of *Alternaria radicina* in *Daucus carota* (carrot) seed by blotter method

**Host:** *Daucus carota* L.

**Pathogen(s):** *Alternaria radicina* Meier, Drechsler & E.D.Eddy, syn. *Stemphylium radicinum* (Meier, Drechsler & E.D.Eddy) Neergaard

**Prepared by:** ISTA-PDC Method Validation Subcommittee

**Authors:** Sheppard, J.W., Cockerell, V & Roberts, S.J.

### Revision history

Version 1.0, 2003-01-01

Version 1.1, 2013-01-01: Definition of sample size

Version 1.2, 2014-01-01: Addition of positive control; common name of host added

Version 1.3, 2017-01-01: Reporting results revised

Version 1.4, 2021-01-01: Sample size added and Methods revised

Version 1.5, 2024-01-01: Sample size and Figures revised

Version 1.6, 2025-01-01: Methods and figure legends revised

### Background

This method was originally published in the *ISTA Handbook of Seed Health Testing* in November 1964 as S.3. No. 5 and was revised by Gambogi (1987). It was incorporated into the *Annexe to Chapter 7: Seed Health Testing Methods* as method 7-002 (Sheppard & Cockerell, 2002). It has been renumbered (7-002a) and slightly modified following studies conducted using six seed lots in 11 laboratories by the International Seed Health Initiative for Vegetable Crops, ISF (ISHI-Veg) in 1999 and 2001 (Van Bilsen, 2003). The studies compared blotter and malt agar methods and concluded that the two were equivalent. Note that seeds can be simultaneously tested for the presence of *Alternaria dauci* using the same method (see method 7-001a).

### Treated seed

Seed treatments may affect the performance of this test. It should only be performed on untreated seed.

### Sample size

The sample size (total number of seeds to be tested) depends on intended use, the maximum acceptable infection level and the analytical sensitivity of the method. The minimum sample size should be 400 seeds.

### Materials

**Reference material:** reference cultures or other appropriate material

**Substrate:** blotters or filter papers, 90 mm, circular (e.g. Whatman No 1 or equivalent), free from micro-organisms and inhibitors (3 per plate)

**Plates:** 90 mm sterile Petri dishes, one per ten seeds

**Incubator:** operating at  $20 \pm 2$  °C, equipped with timer-controlled near-ultraviolet lights (NUV, peak at 360 nm, e.g. colour number 08, Philips; BLB Sylvania)

**Freezer:** operating at  $-20 \pm 2$  °C.

### Sample preparation

It is vital to exclude any possibility of cross-contamination between seed samples. This can be achieved by swabbing/spraying equipment and gloved hands with 70 % ethanol.

### Methods

Critical control points are indicated by CCP.

1. Pretreatment: None.
2. Place three 90 mm filter papers in each plate and soak with sterile distilled/deionised water. Drain away excess water.
3. Plating
  - 3.1 Aseptically place 10 seeds, evenly spaced (CCP), on the surface of the filter paper in each plate.
  - 3.2 Positive control (reference material): Aseptically place seeds evenly spaced (CCP) onto the surface of the filter paper in an appropriate number of plates to obtain the reference culture, or plate a reference culture on media. The number of plates required will depend on the level of contamination of the positive-control seed lot.
4. Incubate for 3 d at  $20 \pm 2$  °C in the dark.

5. Transfer plates to freezer and maintain at  $-20 \pm 2$  °C for 24 h.
6. After freezing, incubate for 6 d at  $20 \pm 2$  °C with alternating 12 h periods of darkness and light, preferably NUV (ISTA, 1984; Tempe, 1968). Plates should be approx. 25 cm below the lights and should not be stacked.
7. Examine seeds under a stereoscopic microscope at  $\times 30$  for fungal growth and up to  $\times 80$  for identification of conidia. Compare with positive control. Record the number of infected seeds in each plate. Conidiophores are simple or occasionally branched, arising usually singly from the surface of the seed, on the emerging radicle or on aerial mycelium (Fig. 1). Conidia are produced singly or in groups of 2 or rarely 3, ellipsoidal or barrel shaped, with little evidence of beak, up to 75  $\mu\text{m}$  long, olivaceous brown, (Ellis, 1971). Under the stereoscopic microscope, conidia appear blackish and glossy (Fig. 1).

## General methods

**Checking tolerances:** Tolerances provide a means of assessing whether or not the variation in results within or between tests are sufficiently wide as to raise doubts about the accuracy of the results. Suitable tolerances, which can be applied to most direct seed health tests, can be found in Table 5B Part 1 of Chapter 5 of the ISTA Rules, or Table G1 in Miles (1963).

**Reporting results:** The result of a seed health test should indicate the scientific name of the pathogen detected and the test method used. When reported on an ISTA Certificate, results are entered under 'Other Determinations'.

The report must indicate the number of seeds tested. In the case of a negative result (pathogen not detected), the results must be reported as 'not detected'.

In the case of a positive result, the report must indicate the percentage of infected seeds.

## Quality assurance

### Specific training

This test should only be performed by persons who have been trained in fungal identification or under their direct supervision.

## Critical control points (CCP)


Spreading hyphae may lead to contamination of other seeds. Seeds must therefore be spaced at least 20 mm from each other, i.e. no more than 10 seeds per 90 mm Petri dish (Step 3).

## References

- Ellis, M. B. (1971). *Dematiaceous Hyphomycetes*. Commonwealth Mycological Institute, Kew, Surrey, England. 608 pp.
- Gambogi, P. (1987). *ISTA Handbook of Seed Health Testing, Working Sheet No. 5: Alternaria radicina on Daucus carota*. ISTA, Zurich, Switzerland.
- International Seed Testing Association (1984). *Report of 18th International Seminar on Seed Pathology*, 9–16 July, 1984, Washington State University, Puyallup, U.S.A.
- Miles, S. R. (1963). Handbook of tolerances and of measures of precision for seed testing. *Proceedings of the International Seed Testing Association*, **28** (3), 525–686.
- Roberts, S. J., Phelps, K., Taylor, J. D. & Ridout, M. S. (1993). Design and interpretation of seed health assays. In *Proceedings of the First ISTA Plant Disease Committee Symposium on Seed Health Testing, Ottawa, Canada* (ed. J. W. Sheppard), pp. 115–125, Agriculture Canada, Ottawa, Canada.
- Sheppard, J. W. & Cockerell, V. (2002). Detection of *Alternaria dauci* on *Daucus carota* (Carrot). *International Rules for Seed Testing. Annexe to Chapter 7: Seed Health Testing Methods*. 7-002.
- Tempe, J. de (1968). The quantitative effect of light and moisture on carrot seed infections in blotter medium. *Proceedings of the International Seed Testing Association*, **33**, 547–553.

## Validation references

- ISTA (2003). Report of a comparative test on *Alternaria dauci* and *Alternaria radicina* on carrot seed. *Method Validation Reports*. International Seed Testing Association, Bassersdorf, Switzerland.

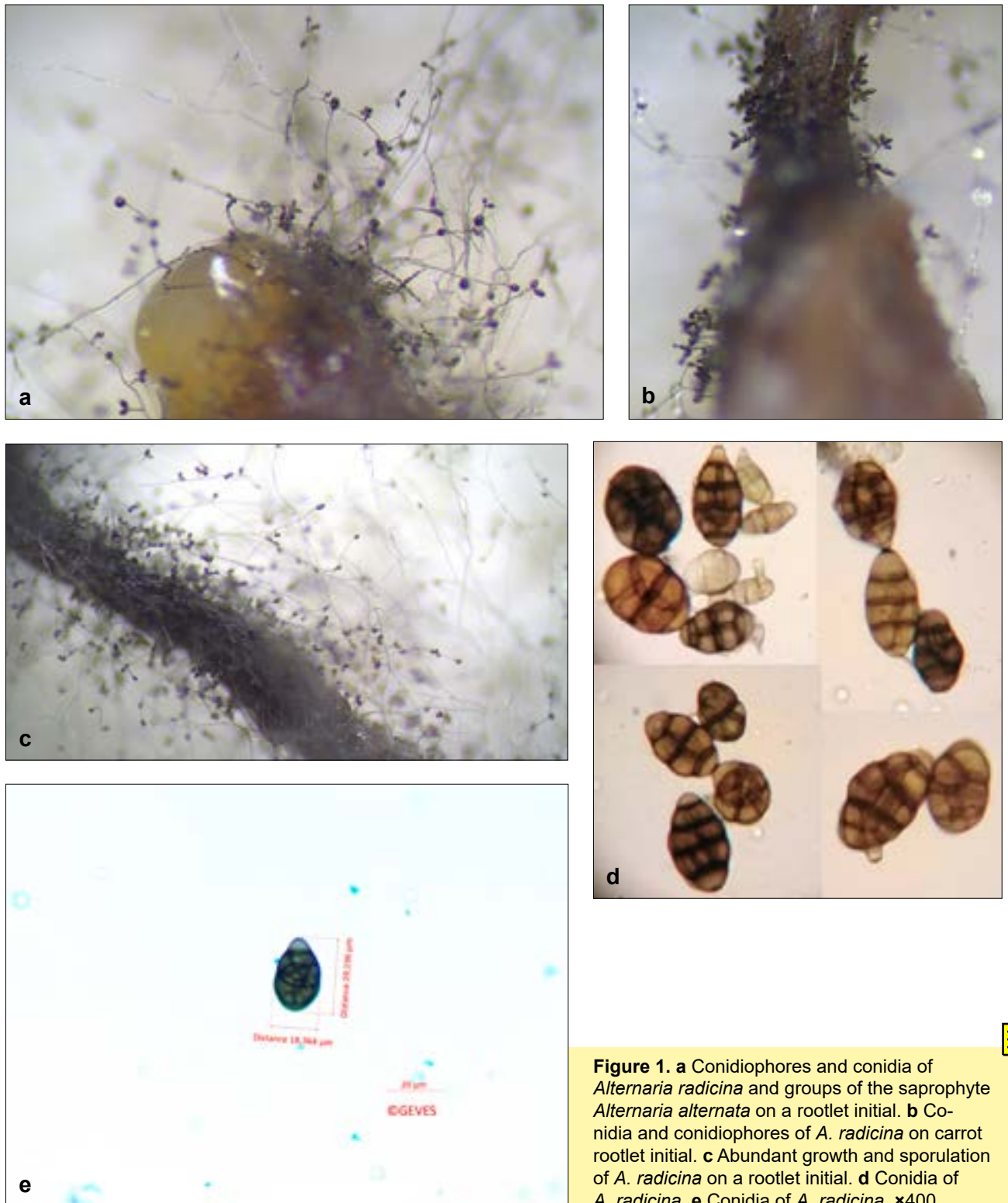

**Figure 1.** a Conidiophores and conidia of *Alternaria radicina* and groups of the saprophyte *Alternaria alternata* on a rootlet initial. b Conidia and conidiophores of *A. radicina* on carrot rootlet initial. c Abundant growth and sporulation of *A. radicina* on a rootlet initial. d Conidia of *A. radicina*. e Conidia of *A. radicina*. ×400.

7-002a: *Alternaria radicina* in *Daucus carota* (carrot) by blotter method

